

OPENLAB VERIFIED · Report identity and source-transcribed quantity. Compound purity and analytical events retain their separate source scopes.

ENDURASHRED

Lab record

Batch	Laboratory	Report	Analysed
OL2212	Janoshik Analytical	29644	16 April 2023

Powered by EvidenceOS, connecting this Olympus Labs UK product and batch to its laboratory results and source evidence.

LGD-4033 · Concentration

8.68 mg

Reported sample; the laboratory gives mg without a per-capsule denominator.

Janoshik report 29644 · OL2212

LGD-4033 · Purity

>99%

Compound raw material

Olympus Labs UK · compound purity statement

MK-2866 · Concentration

7.94 mg

Reported sample; the laboratory gives mg without a per-capsule denominator.

Janoshik report 29644 · OL2212

MK-2866 · Purity

>99%

Compound raw material

Olympus Labs UK · compound purity statement

REPORTED

Compound identity

LGD-4033 · MK-2866

REPORTED

Reported quantity

LGD-4033 8.68 mg · MK-2866 7.94 mg

REPORTED

Compound purity

LGD-4033 >99% · MK-2866 >99%

REPORTED

Laboratory source

Janoshik

REPORTED

Batch connection

OL2212

REPORTED

Sample receipt

2023-04-14

The batch report and supporting compound analyses retain their separate dates and source identities. Compound raw-material purity is a supporting observation; production stage is unspecified.



Compound analysis

Supporting LC/MS analysis · 16 August 2022

LC/MS 16 August 2022 DAD1 B, Sig=220,8 Ref=off

LC/MS Report

File D:\CHEM32\DATA\202208\20220816-KEHU-1 2022-08-16 13-53-55\O-20220816-LCMS.D
Injection Date : 16-Aug-22, 13:54:23 Tgt Mass (EZ) :
Sample Name : Location : P1-F-07
Acq. Operator : Inj : 1
Spec. Reported : MS Integration Inj Volume : 1.5 ul
Acq. Method : D:\CHEM32\DATA\202208\20220816-KEHU-1 2022-08-16 13-53-55\401-CZ-RK.M
Analysis Method : C:\Chem32\1\METHODS\401-GP50-KEHU.M
Sample Info :
Method Info : Instrument : Agilent 1200-6110 MS
Mobile Phase : A=0.1%TFA in water, B=Acetonitrile
Gradient : 5-95% B in 4 min, 95%B in isocratic for 3min
Flow Rate : 0.6 mL/min
Column : XDB-C18 2.1*50mm, 3.5um
Column temperature : 35°C

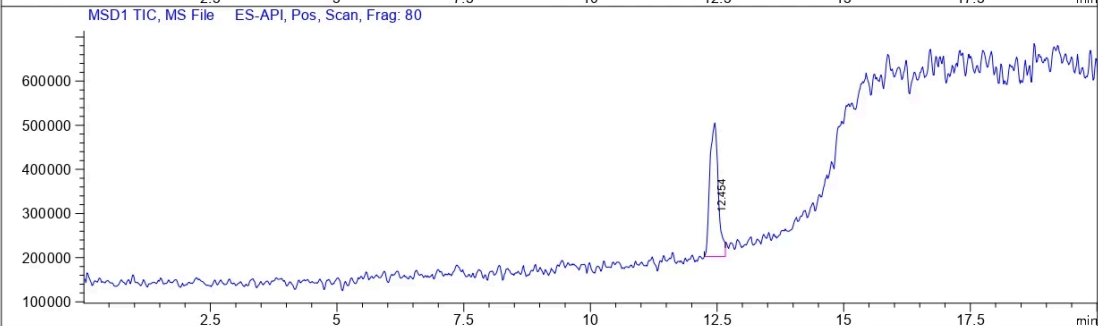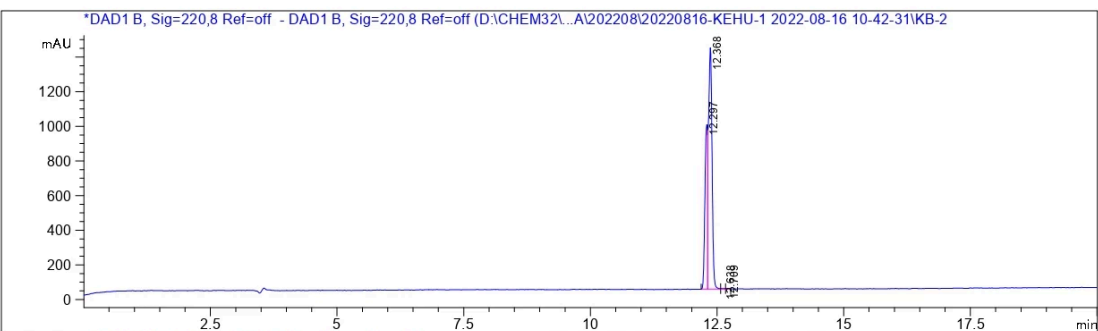

Integration Results for DAD1 B, Sig=220,8 Ref=off

RetTim	Width	Area	Height	Area%
12.30	0.06	3375.57	949.33	33.15
12.37	0.07	6783.10	1392.44	66.62
12.64	0.06	11.87	2.50	0.12
12.71	0.05	11.99	3.48	0.12

Compound analysis


Ostarine LC/MS · companion spectrum · 16 August 2022

LC/MS 16 August 2022 Original source spectrum

LC/MS Report

Ret. Time: 12.45

<<<< POSITIVE SPECTRA >>>>

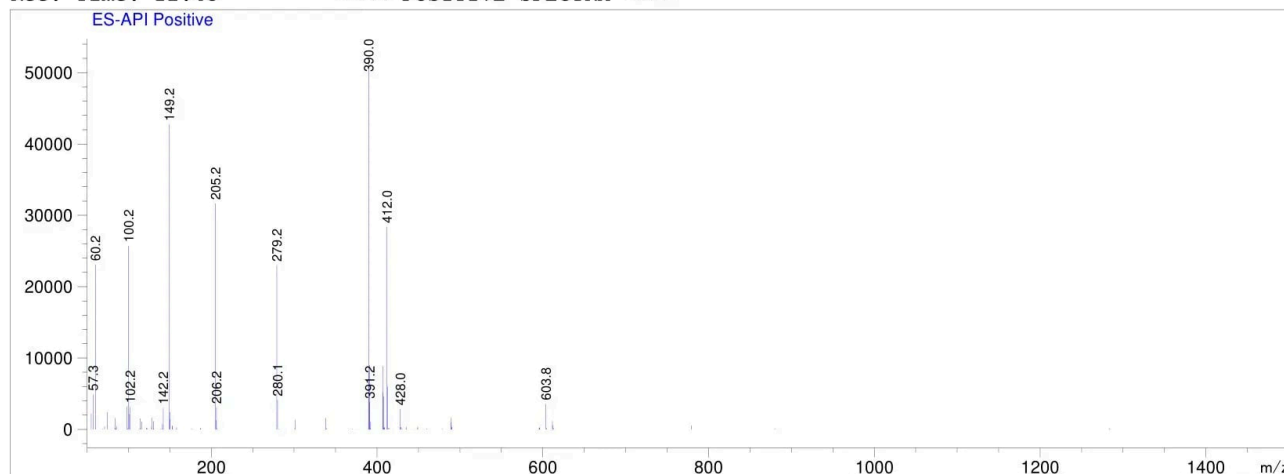

Original analytical page · Supporting compound evidence



Compound analysis

LGD-4033 · LC/MS chromatogram · 2022-08-16

LC/MS 16 August 2022 DAD1 B, 220 nm

LC/MS Report

File D:\CHEM32\DATA\202208\20220816-KEHU-1 2022-08-16 10-42-31\L-20220816-LCMS.D
Injection Date : 16-Aug-22, 11:38:29 Tgt Mass (EZ) :
Sample Name : Location : P1-E-07
Acq. Operator : Inj : 1
Spec. Reported : MS Integration Inj Volume : 3 ul
Acq. Method : D:\CHEM32\DATA\202208\20220816-KEHU-1 2022-08-16 10-42-31\401-CZ-RK.M
Analysis Method : C:\Chem32\1\METHODS\401-GP50-KEHU.M
Sample Info :
Method Info : Instrument : Agilent 1200-6110 MS
Mobile Phase : A=0.1%TFA in water, B=Acetonitrile
Gradient : 5-95% B in 4 min, 95%B in isocratic for 3min
Flow Rate : 0.6 mL/min
Column : XDB-C18 2.1*50mm, 3.5um
Column temperature : 35°C

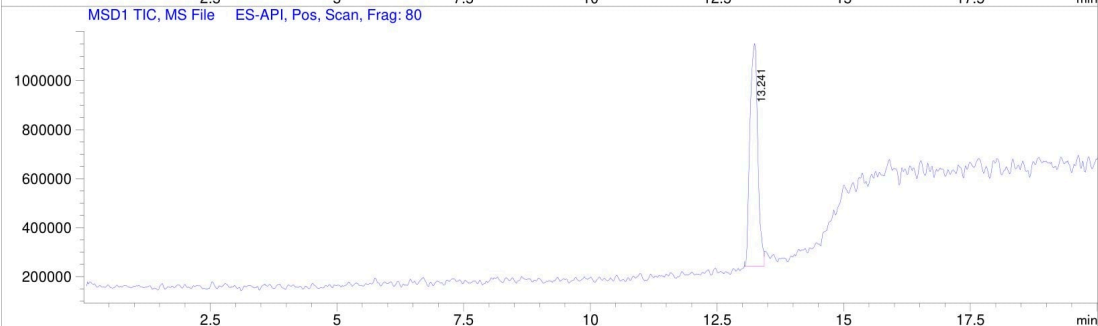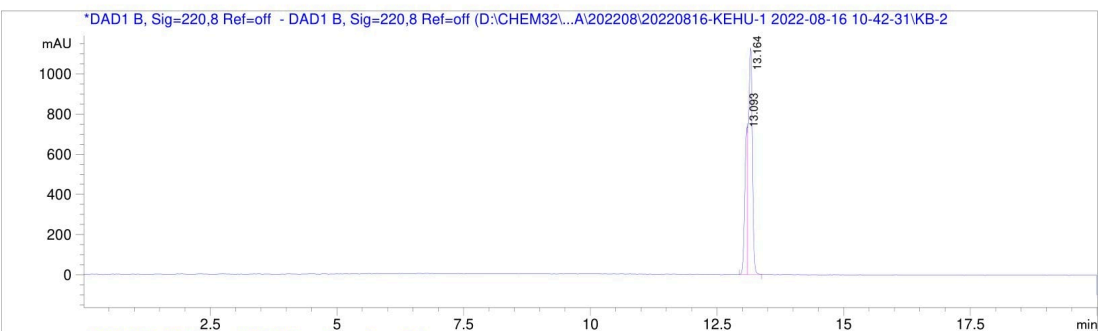

Integration Results for DAD1 B, Sig=220,8 Ref=off


RetTim	Width	Area	Height	Area%
13.09	0.06	2673.60	734.40	30.74
13.16	0.08	6023.74	1128.59	69.26

Compound analysis

LGD-4033 · mass spectrum · 2022-08-16

LC/MS 16 August 2022 Original source spectrum

LC/MS Report

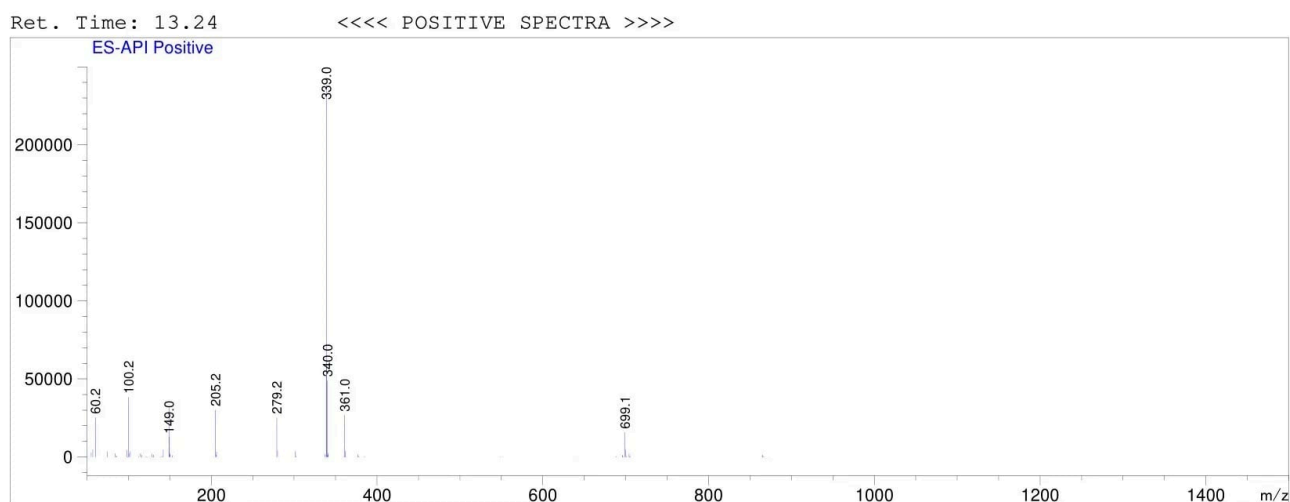

Original analytical page · Supporting compound evidence

Sample analysis

ENDURASHRED · PDA Ch5 227 nm · 2023-04-14

Chromatography 14 April 2023 PDA Ch5 227 nm




Analysis Report

<Sample Information>

Sample Name : 1275-2
Method Filename : Raws - blind common anabolic steroid and aromatase inhibitors screening.lcm
Vial # : 3-2
Injection Volume : 0.5 uL
Date Acquired : 14/04/2023 10:31:04
Date Processed : 14/04/2023 10:35:39
Sample Type : Unknown
Acquired by : Janoshik s.r.o.
Processed by : Janoshik s.r.o.

<Chromatogram>

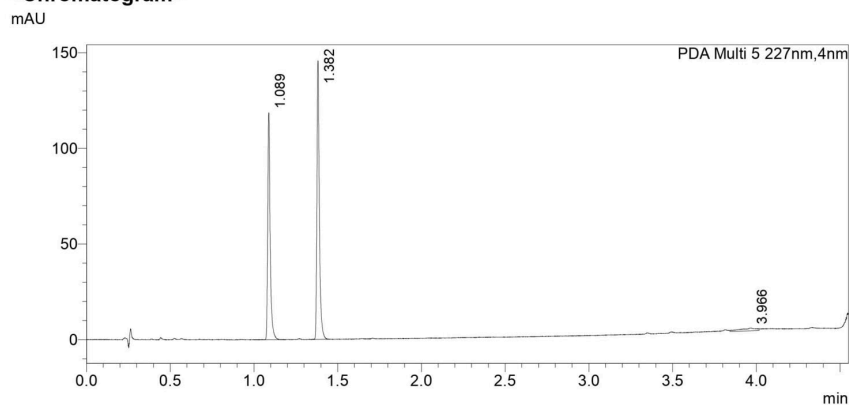

<Peak Table>

PDA Ch5 227nm

Peak#	Ret. Time	Concentration	Unit	Name
1	1.089	7.94	mg	MK-2866
2	1.382	8.68	mg	LGD-4033
3	3.966	0.00		
Total		16.62		

info@janoshik.com

Original analytical page · Report 29644

Report 29644 · OL2212

Sample received 14 April 2023 · Analysed 16 April 2023

TEST REPORT



E-mail: info@janoshuk.com
Web: www.janoshuk.com

Task Number #29644

Testing ordered > 31 MAR '23

Sample received > 14 APR '23

Client support@olympuslabs.uk

Sample LGD-4033 8.5mg / MK-2866 8mg

Manufacturer Olympus Labs UK

Batch OL2212

Sample description >

Pictures not included.

Tests requested >

Qualitative and quantitative analysis of non-AAS sample.

Results >



MK-2886	7.94 mg
LGD-4033	8.68 mg

Comments >

Analysis conducted > 16 APR 2023

Signature >

Verify this test at www.janoshuk.com/verify/ with the following unique key

TST4ESLRXE85

Sources and record edition

Laboratory dates

Sample received: 14 April 2023

Analysed: 16 April 2023

Supporting LC/MS analysis · 16 August 2022

[View source](#)

bc73739c33804ac9bc8d087f9cb6c3ad21b2693bc7f1ef5761fdaf25f1ec72bb

LGD-4033 · LC/MS chromatogram · 2022-08-16

Source image included in this edition.

a7b8c98e5528d64a429e27968ed4ea9f177868e5aa5314f488b074d77d1c19c7

Janoshik report 29644

[View source](#)

22440904d51c4c951fc977dea695d4844a394260e7d20efc9458ac7a2192cb7d

ENDURASHRED · PDA Ch5 227 nm · 2023-04-14

[View source](#)

c379abd2d2749e0820ef9133111102e05f42e3e0a3cf0cd8e9616962d45f3033

Janoshik · 29644

[View source](#)

b8b3705c14b3a22fabd187a5ad246b7c232388ddca92ea50c2081abb6b647531

Ostarine LC/MS · companion spectrum · 16 August 2022

Source image included in this edition.

08f313a408daeb5bf8a95d841bdadb2043bbe6c39c126c46cb4d82cbcf272645

LGD-4033 · mass spectrum · 2022-08-16

Source image included in this edition.

18e88391bb792a61207e3c39a50a0fd4fa3637eba556eac40467308cd8454b85

ENDURASHRED · publisher purity statement

[View source](#)

c90501906a946ca8b065f8b90ea794a3a720568fa48647f865a965a24b2c8329

Edition

openlab.v2.f270d62e05b398f1

Record digest:

47fc5474ae732bc68f2800f0dcc8b21b0ea68c36a7ef69d5606bb44dff72bde5